

This Page Is Inserted by IFW Operations  
and is not a part of the Official Record

## **BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning documents *will not* correct images,  
please do not report the images to the  
Image Problem Mailbox.**

DATE: \_\_\_\_\_ PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_

SUBJECT Inyte clones from LawrencevilleFor hKCNQ5, see p. 72

Steve Dworetzky  
Neuroscience  
Bristol-Myers Squibb  
5 Research Parkway  
Wallingford, CT 06450

Steve Dworetzky,

In p10cy vector

Enclosed you will find a LB + Amp petri dish streaked with the  
1858095 and 3458089. These all originate from a single colony and are stored as glycerol stocks. If you  
have any questions, please do not hesitate to contact me.

already haveNaCh - ASIC-like familyhKCNQ4 (partial)(partial)  
hKCNQ5

chloride ch.  
clones 3447387, 4970006, 3532089  
chloride ch.

Sincerely,

*Tiziano DiPaolo*

Tiziano DiPaolo  
Applied Genomics  
Core Sequencing Facility  
Bristol-Myers Squibb  
Lawrenceville, NJ  
609-252-3991  
dipaolot@bms.com

30 These clones have already been sequenced -

1858095 - we already had - NaCh - ASIC-like family

3447387, 3532089 - are chloride channels - give to Mark Thompson  
to work on -

35 4970006 - is a partial clone: hKCNQ5 } potassium channels  
3458089 - is a partial clone: hKCNQ4 } Curled

SIGNED

DATE

WITNESSED AND  
UNDERSTOOD BY:

DATE

CROSS REFERENCES:

SUBJECT WCC NQ5

See page 42 for Inayate clone info

EST - from Lawrenceville - 4790006 - partial -

hkcna5

Date:

From: "Steven I Dworetzky" <dworetzs@bms.com> Internal  
Organization: Bristol-Myers Squibb  
To: "Trojnacki, Joanne T" <trojnacj@bms.com>

Clone: 4970006  
Library: kidney epithelial transf embryo line, 293-EBNA, lg cDNA  
Putative Annotation: KCNQ Family Member ??  
Vector: pINCY  
Sequence:

sequence:

CTTTTTTTTTTTTTTTTTTTTTTTGGAGATGGAGTCTCTCTCTCTCACCGAGACTTGCC  
CTAAATAAATTTATGGCACAATACAGTAGTGAAGTAAACCAAGCTCAAAAGAGAAATCTGCTGT  
CTCAACAATGACCGCTCAGTACACTAGTACTAGTAACTCTAGGCATCAATCATATATTAC  
TCATTAAATCTGCACAACTACTCTTAGAGACTGGTATTACTAATATTCCCCCATTTAAAA  
AATGGAGAAATTTGTTTATTCTTACGACCAATGGATAAAATGTTTGAATGAATCTGAAAT  
GTAGTAGTCTGACTCCAATGTCCACATGCCTTTAAACCCTACATGTTTGTGTAGCATCT  
GAAAGATAATTCAGGATACTGTAATGAAAGAGTATGTTTCCACAGCAATTTTGTGAAG  
AAAATAAATTTTGTGCATCTTTTCCACACTGTTCTTCTGGATTTCTTCCCACTTTCCCAT  
CTTATCTCAAATCGAGTTAAAAATATACATACATACATATATACATATGTGTATAGTA  
GGAATGTCATCTATTATAAACCATAGGCCAAAGCATGGGATAAGGTTTTATCCAAGT  
TTTTATTGTGTCATCCCTGGGAACGAGGATTATGATTTGGGCTTCTTTATCTCTGTGTGT  
TCGATACATCACTAGGACAGCAAAATAAGCATATAATACTTTGGGAAGGCCAGGAAAA  
CACTCTGTGCTTGTCTTTTAAACCAGGACAGCAAAACAAACCTATTCTGCTGAAAA  
TTTTCTGGCCCAACATGAAGTGTGAAGAACCTTGACTGATCATCCGACTCATTTAA  
ATCTATTGAATGTATTCTAGTACTCAATGTTTCATATAAGTAATAATAATCATTTG  
GGAATATATCTGACCTTACCTCAGGACAGCAAAATGATTCATGTGCTGCTAGTATT  
GTGAATCACTCTTTTAGGCTCATAGCTTTTGTAGTCACTGTAGGGCTATTACTAACCA  
AAGAGATTTGGTCCCGTTCAGAAAGTAGGTGTCTTTAAACCTAACTGAGCAACAGAA  
TCATCTTATATCTATCCCATCTCAAAAATAATAAATCCAATCTCCCATCTGAGACA  
AAGGCTTTAAGGTGATCTTGTAGCTTTGATCTCTCAATTCACTAGGAGGCTGGGGAG  
AATTAGATGTTTTCTCTGGCTAATACCATTTACTTTCTACAAATAAATCTAGACATTTA  
ATGTATGTTAATGAGATACTGTGTAGTACTTAAACAGAGTAAATGTTCTCAGAAATG  
TGAGTCTCTGCCCTCATCAGCAAAAGTGGCTTATAATACAGAGAAATCTGAAGTCTCT  
TGAGCTATACTGAAATTCCTTAGGCTCAATCTAAGCACTTTGAAGGCCAGGCTCAT  
GGCTCACACCTATAATCCCGCTCTTTGGGAGCCCAAGGCCAGGAGGACTGCTTGAGTCCA  
GGGTTTGAAGACCGGCTGAGCAACATCTGAGACCCCTGTCTTACATATTTTAAAA  
TTAGCTTGGCACTGGTGGCATTACCTGTAATCTAGTCTCTCAAGAGCTGAGGAGGAG  
GATCTGTGAAGTGTGTGGGGCTGTGGCTTGGGTGGCGGTATTATTGTTTGTGAATGGT  
GTTGGTGGCTGGCCACTGTGAGCCATGCCTTGACTAATTGGGCACTGTGTTGCTGACT  
GTGACTGTAGTGGGCTAGCGGCTGGAAGAGTCTGGGCATGAGTCACTTTGGCGTCAGAA  
GCACTGACGGCTCTCGAGTTTGTGGCACTAGTTGACTGGATGAGGACGACTGTTTTC  
TGGCGAAGCCCGAAGATCTTTGCTATCCACAGGCGCTTTGATAGTCAGATGTGTTTCACA  
TCCAAAGGTGGGATCTTGAATGAAGCCAAAGCCAGGCGCTGGGACGAGGCTTTCCGAAG  
GACCTCTGATAGATCTCTAGTAGGCAGTCTGAGCTGGATTTCTATGGACTGTACCTGTTT  
TCAACCTTGACCCACGCGGAGCATGACTGAGATGCTGTGGTCTCATGTTCTGCTGCT  
TATTTTCTCTCGGCTCTTCTATCTGATGTGATTTGCCCTTTTCCCTGGAGCGGAAT

sense paper

214846.5

5' RACE primer

2157 5 22 1

5' KACG primer: 5' GC: TAA CAT CAG ACC ACA G 3' Tm = 58°  
6521

5' RAG prim. : 5' CTG GGT CCG GTG GTC AAG GTT 3' Tm = 68°  
65°

5' AACGTTT 3' 5' GGA CTG CTT ACT TGA CAT CTA TC 3' Tm: 68  
65°

h. kDNA 5' sense primer: 5' CTGGATAG CAG CCA CTG TTT 3' Tm: 62

1 of 1

3/9

: Design primers to use for 5' RAGE: (see GIBCO's 5' RAGE protocol)  
Also one primer to use in PCR first - to see if hKCNAS is actually present in HEK 293 cells.

35

**SIGNED**

DATE \_\_\_\_\_

WITNESSED AND  
UNDERSTOOD BY:

DATE \_\_\_\_\_

CROSS REFERENCES:



J\_05\_42 + CTTCTAATAGCAGTGTGTTTCGCAAAACTCAGGGTAATAATTTTTGCCACGCTGCACCT  
consensus GCTTCAAATAGCAGTGTGTTTCGCAAAACTCAGGGTAATAATTTTTGCCACGCTGCACCT

J\_05\_42 + AGAAGCTCCGTTTCCTACAGATCTCCCGCATGGTGCGCATGGACGAAGGGGAGGCACT  
consensus AGAAGCTCCGTTTCCTACAGATCTCCCGCATGGTGCGCATGGACGAAGGGGAGGCACT

J\_05\_42 + TGGAAATTACTGGGTCAGTGGTTATGTCTCACAGAAGGAATTAATCACAGCTTGGTAC  
consensus TGGAAATTACTGGGTCAGTGGTTATGTCTCACAGAAGGAATTAATCACAGCTTGGTAC

J\_05\_42 + ATAGGATTTTGGTCTTAATTTTTTCGCTTTTCCTGTCTATCTGCTGGGAAAAGGATGCC  
J\_05\_67 + KGGGTGGAAGAGGATGCG  
consensus ATAGGATTTTGGTCTTAATTTTTTCGCTTTTCCTGTCTATCTGCTGGGAAAAGGATGCC

J\_05\_42 + A-ATAAAGATTTTCTACATATGCAGATGCTCTCTGGGGGACAATATCAATGCAAC  
J\_05\_67 + AATAAAGATTTTCTACATATGCAGATGCTCTCTGGGGGACAATATCAATGCAAC  
consensus AATAAAGATTTTCTACATATGCAGATGCTCTCTGGGGGACAATATCAATGCAAC

J\_05\_42 + TATTGGCTATGGAGCAAAACTCCCCAACTTGGCTGGGAAGATGCTTTCTGCAGGCTT  
J\_05\_67 + TATTGGCTATGGAGCAAAACTCCCCAACTTGGCTGGGAAGATGCTTTCTGCAGGCTT  
J\_05\_61 + .BGG-GATTGCTTTCTGCGRNCTT  
consensus TATTGGCTATGGAGCAAAACTCCCCAACTTGGCTGGGAAGATGCTTTCTGCAGGCTT

J\_05\_42 + TGCACTCTCTGGCAATTTCTTTCTTTCGACTTCTCTGGGGGATCTTGGCTCAGGTTTTCG  
J\_05\_67 + TGCACTCTCTGGCAATTTCTTTCTTTCGACTTCTCTGGGGGATCTTGGCTCAGGTTTTCG  
J\_05\_61 + TGCACTCTCTGGCAATTTCTTTCTTTCGACTTCTCTGGGGGATCTTGGCTCAGGTTTTCG  
consensus TGCACTCTCTGGCAATTTCTTTCTTTCGACTTCTCTGGGGGATCTTGGCTCAGGTTTTCG

J\_05\_42 + ATTTAAAGATACAGAACACACCGCCAGAAACACTTTGAGAAAGAGGAAGACCAGCTGC  
J\_05\_67 + ATTTAAAGATACAGAACACACCGCCAGAAACACTTTGAGAAAGAGGAAGACCAGCTGC  
J\_05\_61 + ATTTAAAGATACAGAACACACCGCCAGAAACACTTTGAGAAAGAGGAAGACCAGCTGC  
consensus ATTTAAAGATACAGAACACACCGCCAGAAACACTTTGAGAAAGAGGAAGACCAGCTGC

J\_05\_42 + CAACCTCATTCAGTGTGTTTGGCGTAGTTACGCAGCTGATGAGAAATCTGTTTCATTGC  
J\_05\_67 + CAACCTCATTCAGTGTGTTTGGCGTAGTTACGCAGCTGATGAGAAATCTGTTTCATTGC  
J\_05\_61 + CAACCTCATTCAGTGTGTTTGGCGTAGTTACGCAGCTGATGAGAAATCTGTTTCATTGC  
consensus CAACCTCATTCAGTGTGTTTGGCGTAGTTACGCAGCTGATGAGAAATCTGTTTCATTGC

J\_05\_42 + AACCTGGAAGCCACACTTGAAGGCTTGGACACCTCGACGCTTACCAAGAAAGAACAGG  
J\_05\_67 + AACCTGGAAGCCACACTTGAAGGCTTGGACACCTCGACGCTTACCAAGAAAGAACAGG  
J\_05\_61 + AACCTGGAAGCCACACTTGAAGGCTTGGACACCTCGACGCTTACCAAGAAAGAACAGG  
consensus AACCTGGAAGCCACACTTGAAGGCTTGGACACCTCGACGCTTACCAAGAAAGAACAGG

J\_05\_42 + GGAAGCATCAAGCACTTGGGCAAGTGGACATTTGGTCAGAACTAAGTTTAAAGGAGG  
J\_05\_67 + GGAAGCATCAAGCACT-----CAGAACTAAGTTTAAAGGAGG  
J\_05\_61 + GGAAGCATCAAGCACT-----CAGAACTAAGTTTAAAGGAGG  
consensus GGAAGCATCAAGCACT-----CAGAACTAAGTTTAAAGGAGG

J\_05\_42 + AGTGCCATGCGGTAGCCCCAGGGGCCAGAGTATTAAAGGCCGACAGGCTCAGTAGGTA  
J\_05\_67 + AGTGCCATGCGGTAGCCCCAGGGGCCAGAGTATTAAAGGCCGACAGGCTCAGTAGGTA  
J\_05\_61 + AGTGCCATGCGGTAGCCCCAGGGGCCAGAGTATTAAAGGCCGACAGGCTCAGTAGGTA  
consensus AGTGCCATGCGGTAGCCCCAGGGGCCAGAGTATTAAAGGCCGACAGGCTCAGTAGGTA

J\_05\_42 + CAGGAGGTCCTCCCAAGCAGCTGACATACAGGCTGAGGGCAGTCCCACAAAGTGCAGAAAG  
J\_05\_67 + CAGGAGGTCCTCCCAAGCAGCTGACATACAGGCTGAGGGCAGTCCCACAAAGTGCAGAAAG  
J\_05\_61 + CAGGAGGTCCTCCCAAGCAGCTGACATACAGGCTGAGGGCAGTCCCACAAAGTGCAGAAAG  
consensus CAGGAGGTCCTCCCAAGCAGCTGACATACAGGCTGAGGGCAGTCCCACAAAGTGCAGAAAG

J\_05\_42 + CTGAGGCTTCAACGACCGAAGCCGCTTCGGGCGCTTCGTCGGGCTCAAAAGTTCTCAGCC  
J\_05\_67 + CTGAGGCTTCAACGACCGAAGCCGCTTCGGGCGCTTCGTCGGGCTCAAAAGTTCTCAGCC  
J\_05\_61 + CTGAGGCTTCAACGACCGAAGCCGCTTCGGGCGCTTCGTCGGGCTCAAAAGTTCTCAGCC  
consensus CTGAGGCTTCAACGACCGAAGCCGCTTCGGGCGCTTCGTCGGGCTCAAAAGTTCTCAGCC

J\_05\_42 + AAAACCAAGTATAGATGCTGACAGACGCCCCCTGGCACTGATGATGATATATGAGAAAAAGG  
J\_05\_67 + AAAACCAAGTATAGATGCTGACAGACGCCCCCTGGCACTGATGATGATATATGAGAAAAAGG  
J\_05\_61 + AAAACCAAGTATAGATGCTGACAGACGCCCCCTGGCACTGATGATGATATATGAGAAAAAGG  
consensus AAAACCAAGTATAGATGCTGACAGACGCCCCCTGGCACTGATGATGATATATGAGAAAAAGG

J\_05\_42 + ATGCCAGTGTGATGATGATCAGTGGAGAGCTCACCCGCCACTTAAACCTGTCTATTCGAGC  
J\_05\_67 + ATGCCAGTGTGATGATGATCAGTGGAGAGCTCACCCGCCACTTAAACCTGTCTATTCGAGC  
J\_05\_61 + ATGCCAGTGTGATGATGATCAGTGGAGAGCTCACCCGCCACTTAAACCTGTCTATTCGAGC  
consensus ATGCCAGTGTGATGATGATCAGTGGAGAGCTCACCCGCCACTTAAACCTGTCTATTCGAGC

J\_05\_42 + TATCAGAAATTATGAAATTTTCATGTTGCAAAACCGGAAGTTTAAAGGAAACATTACGTCATTA  
J\_05\_67 + TATCAGAAATTATGAAATTTTCATGTTGCAAAACCGGAAGTTTAAAGGAAACATTACGTCATTA  
J\_05\_61 + TATCAGAAATTATGAAATTTTCATGTTGCAAAACCGGAAGTTTAAAGGAAACATTACGTCATTA  
consensus TATCAGAAATTATGAAATTTTCATGTTGCAAAACCGGAAGTTTAAAGGAAACATTACGTCATTA

J\_05\_42 + TGATGTAAGAAGTGTCAATGAAACAATATTTCTGGTGTCTCTGGACATCTGTGTGTAGAT  
J\_05\_67 + TGATGTAAGAAGTGTCAATGAAACAATATTTCTGGTGTCTCTGGACATCTGTGTGTAGAT  
J\_05\_61 + TGATGTAAGAAGTGTCAATGAAACAATATTTCTGGTGTCTCTGGACATCTGTGTGTAGAT  
consensus TGATGTAAGAAGTGTCAATGAAACAATATTTCTGGTGTCTCTGGACATCTGTGTGTAGAT

05/16264 (5893)

J\_05\_42 + TAAAGGCTTCAAACAGCGTGTGATCAAAATCTCTGGAAAGGGCAAAATCACATCAGATAA  
J\_05\_67 + TAAAGGCTTCAAACAGCGTGTGATCAAAATCTCTGGAAAGGGCAAAATCACATCAGATAA  
J\_05\_61 + TAAAGGCTTCAAACAGCGTGTGATCAAAATCTCTGGAAAGGGCAAAATCACATCAGATAA  
J\_05\_47 + TAAAGGCTTCAAACAGCGTGTGATCAAAATCTCTGGAAAGGGCAAAATCACATCAGATAA  
consensus TAAAGGCTTCAAACAGCGTGTGATCAAAATCTCTGGAAAGGGCAAAATCACATCAGATAA

J\_05\_42 + GAAGAGCCGAGAGAAAATACAGAGCAGAACATAGACACCAAGACGATCTCAGTATGCTCGG  
J\_05\_67 + GAAGAGCCGAGAGAAAATACAGAGCAGAACATAGACACCAAGACGATCTCAGTATGCTCGG  
J\_05\_61 + GAAGAGCCGAGAGAAAATACAGAGCAGAACATAGACACCAAGACGATCTCAGTATGCTCGG  
J\_05\_47 + GAAGAGCCGAGAGAAAATACAGAGCAGAACATAGACACCAAGACGATCTCAGTATGCTCGG  
consensus GAAGAGCCGAGAGAAAATACAGAGCAGAACATAGACACCAAGACGATCTCAGTATGCTCGG

J\_05\_42 + TCGGGTGGTCAAGGTGAAAACAGGTACATCTCATAGATCCAAAGCTGGAGTGCCTACT  
J\_05\_67 + TCGGGTGGTCAAGGTGAAAACAGGTACATCTCATAGATCCAAAGCTGGAGTGCCTACT  
J\_05\_61 + TCGGGTGGTCAAGGTGAAAACAGGTACATCTCATAGATCCAAAGCTGGAGTGCCTACT  
J\_05\_47 + TCGGGTGGTCAAGGTGAAAACAGGTACATCTCATAGATCCAAAGCTGGAGTGCCTACT  
consensus TCGGGTGGTCAAGGTGAAAACAGGTACATCTCATAGATCCAAAGCTGGAGTGCCTACT

J\_05\_42 + AGACATCTATCAACAGGTCTCTCGGAAAGGCTCTGCCCTCAGCCCTCGCTTGGCTTCAATT  
J\_05\_67 + AGACATCTATCAACAGGTCTCTCGGAAAGGCTCTGCCCTCAGCCCTCGCTTGGCTTCAATT  
J\_05\_61 + AGACATCTATCAACAGGTCTCTCGGAAAGGCTCTGCCCTCAGCCCTCGCTTGGCTTCAATT  
J\_05\_47 + AGACATCTATCAACAGGTCTCTCGGAAAGGCTCTGCCCTCAGCCCTCGCTTGGCTTCAATT  
consensus AGACATCTATCAACAGGTCTCTCGGAAAGGCTCTGCCCTCAGCCCTCGCTTGGCTTCAATT

J\_05\_42 + CCAGATGCCACCTTTTGAATGTGAACAGACATCTGACTATCAAGGCCCTGTGGATAGCAA  
J\_05\_67 + CCAGATGCCACCTTTTGAATGTGAACAGACATCTGACTATCAAGGCCCTGTGGATAGCAA  
J\_05\_61 + CCAGATGCCACCTTTTGAATGTGAACAGACATCTGACTATCAAGGCCCTGTGGATAGCAA  
J\_05\_47 + CCAGATGCCACCTTTTGAATGTGAACAGACATCTGACTATCAAGGCCCTGTGGATAGCAA  
consensus CCAGATGCCACCTTTTGAATGTGAACAGACATCTGACTATCAAGGCCCTGTGGATAGCAA

J\_05\_67 + AGATCTTTTGGGTTCCGCAACAAACAGTGGCTGCTTATCCAGATCAACTAGTGCCAACT  
J\_05\_61 + AGATCTTTTGGGTTCCGCAACAAACAGTGGCTGCTTATCCAGATCAACTAGTGCCAACT  
J\_05\_47 + AGATCTTTTGGGTTCCGCAACAAACAGTGGCTGCTTATCCAGATCAACTAGTGCCAACT  
consensus AGATCTTTTGGGTTCCGCAACAAACAGTGGCTGCTTATCCAGATCAACTAGTGCCAACT

J\_05\_67 + CTCGAGAGGCGTCGAGTTCATTCTGACGCCAAATGAGTTCAGTGCCCAAGCTTTCTACGC  
J\_05\_61 + CTCGAGAGGCGTCGAGTTCATTCTGACGCCAAATGAGTTCAGTGCCCAAGCTTTCTACGC  
J\_05\_47 + CTCGAGAGGCGTCGAGTTCATTCTGACGCCAAATGAGTTCAGTGCCCAAGCTTTCTACGC  
consensus CTCGAGAGGCGTCGAGTTCATTCTGACGCCAAATGAGTTCAGTGCCCAAGCTTTCTACGC

J\_05\_67 + GCTTAGGCCCTACTATGCAAGCTCAAGCAACACAGGTCGCAATAGTCAAAAGCGATGGCTC  
J\_05\_61 + GCTTAGGCCCTACTATGCAAGCTCAAGCAACACAGGTCGCAATAGTCAAAAGCGATGGCTC  
J\_05\_47 + GCTTAGGCCCTACTATGCAAGCTCAAGCAACACAGGTCGCAATAGTCAAAAGCGATGGCTC  
consensus GCTTAGGCCCTACTATGCAAGCTCAAGCAACACAGGTCGCAATAGTCAAAAGCGATGGCTC

J\_05\_67 + AGCACTGGCCAGCCACCAACCAATTCGCAACCAAAATAATACGGCAACCAAGCCAGCAAGC  
J\_05\_61 + AGCACTGGCCAGCCACCAACCAATTCGCAACCAAAATAATACGGCAACCAAGCCAGCAAGC  
J\_05\_47 + AGCACTGGCCAGCCACCAACCAATTCGCAACCAAAATAATACGGCAACCAAGCCAGCAAGC  
consensus AGCACTGGCCAGCCACCAACCAATTCGCAACCAAAATAATACGGCAACCAAGCCAGCAAGC

J\_05\_67 + CCCAACAACTTTACAGATCCCACTCTCTCTCCAGGCATCAAGCATCTGCCAGGCCAGA  
J\_05\_61 + CCCAACAACTTTACAGATCCCACTCTCTCTCCAGGCATCAAGCATCTGCCAGGCCAGA  
J\_05\_47 + CCCAACAACTTTACAGATCCCACTCTCTCTCCAGGCATCAAGCATCTGCCAGGCCAGA  
consensus CCCAACAACTTTACAGATCCCACTCTCTCTCCAGGCATCAAGCATCTGCCAGGCCAGA

J\_05\_67 + AACTCTGCACCTTAACCTTGAGGCTTACAGGAAAGCAATTTCTGAGCTCAGCACTGGCT  
J\_05\_61 + AACTCTGCACCTTAACCTTGAGGCTTACAGGAAAGCAATTTCTGAGCTCAGCACTGGCT  
J\_05\_47 + AACTCTGCACCTTAACCTTGAGGCTTACAGGAAAGCAATTTCTGAGCTCAGCACTGGCT  
consensus AACTCTGCACCTTAACCTTGAGGCTTACAGGAAAGCAATTTCTGAGCTCAGCACTGGCT

J\_05\_67 + TGTGCTCCCAAGGAAAATGTTTCAGGTTGCAAGTCAAAATCTCACCAGGACCGTTCTAT  
J\_05\_61 + TGTGCTCCCAAGGAAAATGTTTCAGGTTGCAAGTCAAAATCTCACCAGGACCGTTCTAT  
J\_05\_47 + TGTGCTCCCAAGGAAAATGTTTCAGGTTGCAAGTCAAAATCTCACCAGGACCGTTCTAT  
consensus TGTGCTCCCAAGGAAAATGTTTCAGGTTGCAAGTCAAAATCTCACCAGGACCGTTCTAT

J\_05\_67 + GAGGAAAAGCTTTGACATGGGAGGAGAACTCTGTTGTCTGTCTGCTCCCATGCTGCCGAA  
J\_05\_61 + GAGGAAAAGCTTTGACATGGGAGGAGAACTCTGTTGTCTGTCTGCTCCCATGCTGCCGAA  
J\_05\_47 + GAGGAAAAGCTTTGACATGGGAGGAGAACTCTGTTGTCTGTCTGCTCCCATGCTGCCGAA  
consensus GAGGAAAAGCTTTGACATGGGAGGAGAACTCTGTTGTCTGTCTGCTCCCATGCTGCCGAA

J\_05\_67 + GGACTTGGGCAAACTTTTGTCTGTGCAAAACCTGATCAGGTGCAGCCAGGAACTGAAATAT  
J\_05\_61 + GGACTTGGGCAAACTTTTGTCTGTGCAAAACCTGATCAGGTGCAGCCAGGAACTGAAATAT  
J\_05\_47 + GGACTTGGGCAAACTTTTGTCTGTGCAAAACCTGATCAGGTGCAGCCAGGAACTGAAATAT  
consensus GGACTTGGGCAAACTTTTGTCTGTGCAAAACCTGATCAGGTGCAGCCAGGAACTGAAATAT

J\_05\_67 + ACAACTTTTCAGGAGTGAAGTGAAGTGGCTCAGAGGCTCCCAAGATTTTACCCCAAGT  
J\_05\_61 + ACAACTTTTCAGGAGTGAAGTGAAGTGGCTCAGAGGCTCCCAAGATTTTACCCCAAGT  
J\_05\_47 + ACAACTTTTCAGGAGTGAAGTGAAGTGGCTCAGAGGCTCCCAAGATTTTACCCCAAGT  
consensus ACAACTTTTCAGGAGTGAAGTGAAGTGGCTCAGAGGCTCCCAAGATTTTACCCCAAGT

J\_05\_67 + GAGGGAATCCAAATTTGTTTAACTGATGAAGAGGTGGGTCGGAAGGACGACGACGACG  
J\_05\_61 + GAGGGAATCCAAATTTGTTTAACTGATGAAGAGGTGGGTCGGAAGGACGACGACGACG  
J\_05\_47 + GAGGGAATCCAAATTTGTTTAACTGATGAAGAGGTGGGTCGGAAGGACGACGACGACG  
consensus GAGGGAATCCAAATTTGTTTAACTGATGAAGAGGTGGGTCGGAAGGACGACGACGACG

$[Q] = 0.5 \text{ probe} = 26 \text{ unit} = 86 \text{ h molar acid}$

conf

DO NOT WRITE IN THIS MARGIN

CROSS REFERENCES:

WITNESSED AND  
UNDERSTOOD BY

SUBJECT *h k c n 25*

→ contd

CROSS REFERENCES:

DATE: \_\_\_\_\_ PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_

SUBJECT hKCRDS

Cont'd from p. 189

Templates for PCR: Use Human Brain cDNA (4x Clontech Marathon-Ready) with the start/stop and start BsrV/stop XbaI primers.  
 → Tube (D) → Tube (E)

5

Also, use OS Plasmid DNA (#42, #64) for: Template

Tube (A)	start RV / (2) antisense	#42 plasmid DNA
Tube (B)	stop / (3) sense	#64 "
Tube (C)	stop XbaI / (3) sense	#64 "

<sup>10</sup> Dilute plasmid DNA (in pCDNA 3.1 + Neo) to 10 ng/μl - use the -  
 #42 - use 2/11 Maxiprep at 2.82 μg/μl  
 #64 - use 2/18 Maxiprep at 2.11 μg/μl

Plasmid PCR:


#42 PCR Tubes (A) - 1140 bp to be amplified (at 1 min / kb) -  
 extend at 68° for 1 min, 30 seconds

#64 PCR tubes (B, C) - 1720 bp to be amplified (at 1 min / kb) -  
 extend at 68° for 2 minutes.

Brain cDNAHuman Brain cDNA PCR:

<sup>20</sup> Tubes (D, E) - Start/stop primers: extend at 68° for 3 minutes.

94° 3 minutes      94° 30 sec      55° 30 sec      68° as indicated above      + final 68° 10 minutes, 25 cycles
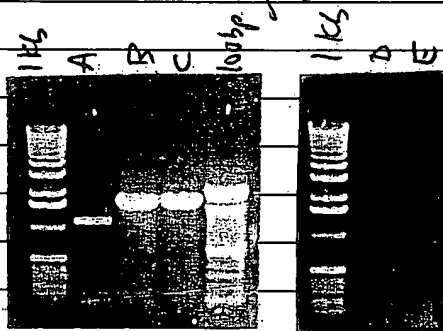
<sup>25</sup> Run 5% agarose gel:

A - expect 1140 bp

B - expect 1720 bp

C - expect 1720 bp

D, E - expect 2200 bp

30

35

SIGNED

DATE

WITNESSED AND  
UNDERSTOOD BY:

DATE

Cont'd →

Barbara G. Gierke

Joe G. Gierke

DO NOT WRITE IN THIS MARGIN

DATE: \_\_\_\_\_ PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_

SUBJECT hkcNOS

→ cont'd

For the brain full length: Tubes D = Kozak start ECDRV / stop  
 E = Kozak start ECDRV / stop XbaI

Have 2700 bp band - shotgun clone 3x into TA:

5 2x vector, 1x ligase, 1x ligase, 3x lvs. Also do TA self-lig, 14° o.m.

(NOTE - will need to prepare pCDNA 3.1(+)-Neo 5' RV/3' XbaI)


Transform 2x, plate out 10x into 50 µg/ml c+amp plate, 37°C.

1 Blue #white

Set up miniprep:

	Blue	White	
10 TA-self-lig	70	5	
start/stop	200	25	#101-110
start/stop XbaI	200	25	#111-120

Wizard Plus Miniprep

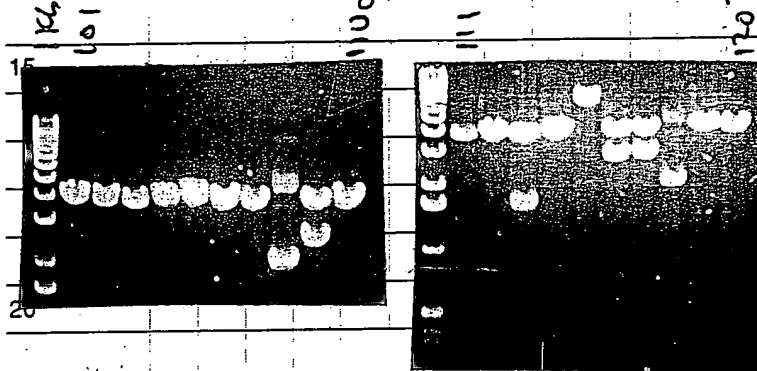

OD, 1:50 :

# 260.0

S# 26

108	0.16884	20 ng/µl
109	0.16934	20 ng/µl
113	0.18014	50 ng/µl
116	0.18264	55 ng/µl
117	0.18304	55 ng/µl
115	-	Assume 400 ng/µl

Sequence:

Req'n #: 18914

CANCEL

Project: HB

For: trojnacj (Joanne Trojnacki)

Owner: trojn

Status: Submitted

Comments: Walli

#109 - M13r - looks good - see start primer, Q5 sequence.  
 - T<sub>7</sub> - weird peaks, but can see stop primer and Q5 sequence.

8 reaction(s) to submit

No.	Template	Primer	Submitter's Comments	Operator's Comments
1.	HB_Q5_contig_PCR_no_109	M13rev		
2.	HB_Q5_contig_PCR_no_109	T7		
3.	HB_Q5_contig_PCR_no_115	M13rev		
4.	HB_Q5_contig_PCR_no_115	T7		
5.	HB_Q5_contig_PCR_no_116	M13rev		
6.	HB_Q5_contig_PCR_no_116	T7		
7.	HB_Q5_contig_PCR_no_117	M13rev		
8.	HB_Q5_contig_PCR_no_117	T7		

See sequencing results on p. 198

#116 - M13r: see stop primer w/ XbaI site, and Q5 sequence.

T<sub>7</sub> - see start primer, RV site, Q5 seq.

#116 is in the wrong orientation in TA, but will cut out Q5 w/ RV-built into start primer - and XbaI built into stop

primer and ligate. Cont'd →

SIGNED

DATE

WITNESSED AND UNDERSTOOD BY:

CROSS REFERENCES:

DO NOT WRITE IN THIS MARGIN



DATE: \_\_\_\_\_

PROJ. NO. \_\_\_\_\_

EXPT. NO. \_\_\_\_\_

SUBJECT: HKCND5

→ cont'd

(See page 192)

Also, have correct sized products for Tubes (A) - start / (2) - #42

1720bp (B) (3) / stop - #64

1720bp (C) (3) / stop Xba - #64

Plasmid  
Template

Plasmid PCR -

5' ~~XXXXXXXXXX~~Run Through Clontech's Nucleotrap PCR Purification - kit - Elute with 40 µl H<sub>2</sub>O,  
for tube (A), 50x for tubes (B), (C); OD, 1:50;

	S#	260.0	280.0	λA/λB	S#	260.0	λ
10 In 40x H <sub>2</sub> O	1 A	0.0028	0.0021	1.294	28	0.0028	1.294
50x	2 B	0.0161	0.0097	1.666	29	0.0161	1.666
50x	3 C	0.0101	0.0056	1.804	30	0.0101	1.804

Use in PCR: Templates A/B, A/C - with start / stop  
dil. each 2 µl - use 2.5 µl or start / stop Xba primers

Use Plat. Tag Hi-Fi

15 94°, 3 min.


94° 30" 55° 30" 68° 10' final  
25 cycles~~XXXXXXXXXX~~

Run 5x on gel:

AB AC B

Looks good.

	260.0	280.0	λA/λB	S#
20	0.0224	0.0178	1.259	56 y/l

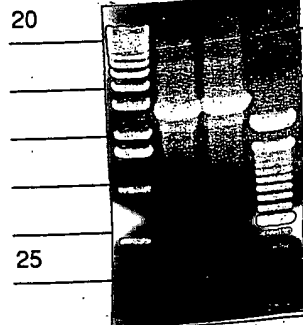

Ligate 2 µl of each into TA vector.

Rest: clean up A/B w/ Nucleotrap PCR Purif. kit.

For tube A/C - clean up w/ PCR purif. kit. Then

Digest w/ EcoRV / XbaI: eluted into 50 µl H<sub>2</sub>O -

save a few µl in case I need it - Rest:

DNA 10x BSA 10x H<sub>2</sub>O Enz.

43 µl 6 6 1 2 µl 2 XbaI, 37°, 2 hr.

Run on gel:

Cut out - clean up w/ Clontech's

Nucleotrap kit - ligate 2 µl into

vector, 140 over weekend -

(Actually was at 160)

From 2/12/95 - 5' RV/3' XbaI at 30 µg/µl

36: Transform 2 µl into DH5α competent cells, plate out 25 µl, 250 µl onto

100 µg/ml Camp plates, 37° env.

35(5) Self-lig - 8 colonies

Expt'l - ~ colonies (250 µl = 30 colonies) pick for minis - #121-140 Cont'd →

Gatop. 196 for minis

SIGNED

DATE

WITNESSED AND  
UNDERSTOOD BY:

DATE

REFERENCES:


DATE: \_\_\_\_\_ PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_

SUBJECT WKEN25

→ cont'd

Need to prepare some pCDNA3.1(+) Neo: 5' EcoRV / 3' Xba I —

Have DNA from 2/26/99 at 1.19 µg/ml — Digest 5 µg:

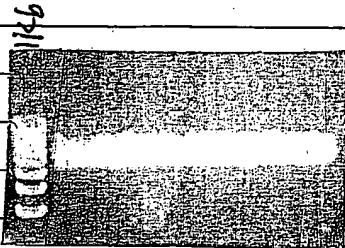DNA <sup>10x</sup> NEB(2) <sup>10x</sup> BSA <sup>10x</sup> H<sub>2</sub>O Eng5 41 5 5 32 2 RV, 2 Xba I, 37°, a few hrs. Run on gel:Clean up Clontech's Nucleotrap kit — Elute 50  
H<sub>2</sub>O, OD, 1:50:

	260.0	280.0	A/A B	S#	260
	0.0228	0.0120	1.899	2.85 µg/µl	
				in 50 µl	
				= 57 ng/µl	


Use 12 in digested w/ 85.

from A/C — RV/Xba I dig'd  
from previous page — use 3 d insert —  
16° o.n.

3/8: Transform 22 → DH5α comp. cells, plate out 25 µl, 250 µl each

100 µg/ml ct amp plates, 37° o.n. Set up miniprep  
 Self-lig (25 µl) — 16 colonies  
 Expt'l (25 µl) — 31 colonies #141-150

Wyzard Plus Mini, digest 22 w/ EcoRV/Xba I:
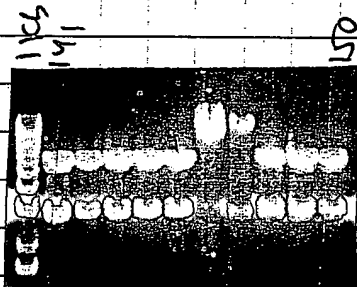
# 260.0  
OD, 1:50:

#141	0.2564	640 ng/µl
142	0.1456	365 ng/µl
143	0.2019	504 ng/µl
146	0.1860	465 ng/µl

Sequence:

8 reaction(s) to submit

No.	Template	Primer	Submitter's Comments	Operator's Comments
1.	Q5_contig_PCR_2700bp_no_141	T7		
2.	Q5_contig_PCR_2700bp_no_141	BGHrev		
3.	Q5_contig_PCR_2700bp_no_142	T7		
4.	Q5_contig_PCR_2700bp_no_142	BGHrev		
5.	Q5_contig_PCR_2700bp_no_143	T7		
6.	Q5_contig_PCR_2700bp_no_143	BGHrev		
7.	Q5_contig_PCR_2700bp_no_146	T7		
8.	Q5_contig_PCR_2700bp_no_146	BGHrev		

See sequencing results in  
Bart 47451, Page 24

Cont'd →

SIGNED

DATE

WITNESSED AND  
UNDERSTOOD BY:

DATE

CROSS REFERENCES:

DATE: \_\_\_\_\_ PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_

SUBJECT hKCNB5


→ cont'd

mini prep #121-140 (see p.194) - Rn A/c (Plasmid PCR) - 5'RV/Xba I dig'd,  
ligated into pCDNA3;

OD, 1:50!

3# 260.0

121	0.1795448ng/λ
132	0.2007500ng/λ
133	0.2202550ng/λ

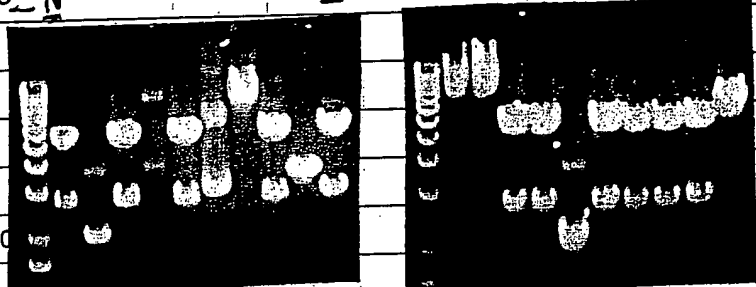

Sequence

Req'n #: 18933 **CANCEL** Project: Q5 Contig Plasmid  
 For: trojnacj (Joanne Trojnacki) Owner: trojnacj (Joanne Troj  
 Status: Submitted Comments: Wallingford

RESULTS:

#132 - Bgt rev: see Xba I site, stop site,  
Q5 sequence.

T7 - don't see start site - BLAST

next page.

6 reaction(s) to submit

No.	Template	Primer	Submitter's Comments	Operator's Comments
1.	Q5_contig_PCR_2700bp_no_121	T7		
2.	Q5_contig_PCR_2700bp_no_121	BGHrev		
3.	Q5_contig_PCR_2700bp_no_132	T7		
4.	Q5_contig_PCR_2700bp_no_132	BGHrev		
5.	Q5_contig_PCR_2700bp_no_133	T7		
6.	Q5_contig_PCR_2700bp_no_133	BGHrev		

#121 - Bgt rev: failed

- T7: data a bit messy, but  
see start primer / start site and  
Q5 sequence.

RESULTS:

Sequencing run # 17931, which contained samples from one or more requests you submitted, is now complete.

Below is the status of the items you requested. You can browse the CSF database for more information at: <http://cyclops.hpw.pri.bms.com:8084/>

\*\*\* PLEASE COME TO THE DNA SEQUENCING LAB NOW TO PICK UP YOUR SAMPLES \*\*\*

Project Status	Requ'n	Item	Template	Primer	Run ID	Lane
Q5_Contig_Plasmid_PCR_2700bp_in_pcDNA3 OK	18933	1	Q5_contig_PCR_2700bp_no_121	T7	17931	12
Q5_Contig_Plasmid_PCR_2700bp_in_pcDNA3 OK	18933	2	Q5_contig_PCR_2700bp_no_121	BGHrev	17931	15
Q5_Contig_Plasmid_PCR_2700bp_in_pcDNA3 OK	18933	3	Q5_contig_PCR_2700bp_no_132	T7	17931	18
Q5_Contig_Plasmid_PCR_2700bp_in_pcDNA3 OK	18933	4	Q5_contig_PCR_2700bp_no_132	BGHrev	17931	21
Q5_Contig_Plasmid_PCR_2700bp_in_pcDNA3 OK	18933	5	Q5_contig_PCR_2700bp_no_133	T7	17931	24
Q5_Contig_Plasmid_PCR_2700bp_in_pcDNA3 OK	18933	6	Q5_contig_PCR_2700bp_no_133	BGHrev	17931	25

\*\*\* End of report \*\*\*

SIGNED

DATE

WITNESSED AND UNDERSTOOD BY:

DATE

Joanne Trojnacki

Ng G/L

DO NOT WRITE IN THIS MARGIN

DATE: \_\_\_\_\_ PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_  
 SUBJECT: hkcns

→ cont'd

B185: #132 (See previous page):

- Sequences producing significant alignments:

	Score (bits)	E Value
ref NM_010611.1  Mus musculus potassium voltage-gated chan...	74	5e-11
ref NM_004518.1  Homo sapiens potassium voltage-gated chan...	74	5e-11
gb AF110020.1 AF110020 Homo sapiens potassium channel (KCNQ...	74	5e-11
gb AF074247.1 AF074247 Homo sapiens neuronal delayed-rectif...	74	5e-11
gb AF033348.1 AF033348 Homo sapiens potassium channel (KCNQ...	74	5e-11
emb Y15065.1 HSCNQ2 Homo sapiens mRNA for voltage gated po...	74	5e-11
dbj D82346.1 D82346 Homo sapiens mRNA for HNSPC, complete cds	74	5e-11
dbj AB000504.1 AB000504 Mus musculus mRNA for alternative s...	74	5e-11

'Q Family' Member

ref|NM\_010611.1| Mus musculus potassium voltage-gated channel, subfamily Q, member 2  
 (Kcnq2), mRNA  
 Length = 2247

Score = 73.8 bits (37), Expect = 5e-11  
 Identities = 55/61 (90%)  
 Strand = Plus / Plus

Query: 51 tgcagaactacctgtacaacgtgctggagagaccccgcggtggcggttcattaccacg 110  
 |||||  
 Sbjct: 316 tgcagaatttcctctacaacgtgctagagcgggcccgcggtggcggttcattaccacg 375

Query: 111 c 111  
 Sbjct: 376 c 376

ref|NM\_004518.1| Homo sapiens potassium voltage-gated channel, KQT-like subfamily,  
 member 2 (KCNQ2) mRNA  
 Length = 7420

Score = 73.8 bits (37), Expect = 5e-11  
 Identities = 55/61 (90%)  
 Strand = Plus / Plus

Query: 51 tgcagaactacctgtacaacgtgctggagagaccccgcggtggcggttcattaccacg 110  
 |||||  
 Sbjct: 272 tgcagaatttcctctacaacgtgctggagcgggcccgcggtggcggttcattaccacg 331

Query: 111 c 111  
 Sbjct: 332 c 332

Score = 50.1 bits (25), Expect = 8e-04  
 Identities = 55/65 (84%)  
 Strand = Plus / Plus

Query: 449 cgcattggaccgaaggggagggcacttggaaattactgggttcagtggtttatgctcacagc 508  
 |||||  
 Sbjct: 670 cgcattggaccggcggggagggcacttggaaagctgctgggtctgtggtctatgccacagc 729

Query: 509 aagga 513  
 |||||  
 Sbjct: 730 aagga 734

Will grow up #121 for maxiprep → Orogen Maxi - elute 250A TE, OD, 1 hr.

260.0 280.0  $\lambda A/\lambda B$  S# 260.00.6073 0.3232 1.879 759  $\mu g$  total

3

= 3.03  $\mu g$  / 1Digest ~ 300 ng w/ EcoRV/Xba, 37°, 1 hr.  
Run on gel -

(Also made glycerol stock)

SIGNED

DATE

WITNESSED AND  
UNDERSTOOD BY:

DATE

CROSS REFERENCES:

DO NOT WRITE IN THIS MARGIN

DATE: \_\_\_\_\_ PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_  
SUBJECT: hKCNDS

→ cont'd

(See p. 193)

I mean time: sequencing results = TA (Hu Brain, Q5 PCR prod (270 bp), #116!

In TA

HB\_Q5\_contig\_PCR\_no\_116\_T7.1791  
[Strand]

EcoRV start

1 GAATTCGGCT TGAATATC CATGAAGGAT GTGGAGTCGG GCCGGGCGAG GGTGCTGCTG AACTCGGCAG CCGCCAGGGG  
GluPheGlyLeu AspIleThr METLysAsp ValGluSerGly ArgGlyArg ValLeuLeu AsnSerAlaAla AlaArgG1  
81 CGACGGCCTG CTACTGCTGG GCACCCGCGC GGCCACGCTT GGTGGCGGCG GCGGTGGCCT GAGGGAGAGC CGCCGGGGCA  
yAspGlyLeu LeuLeuLeuGly ThrArgAla AlaThrLeu GlyGlyGlyGly GlyGlyLeu ArgGluSer ArgArgGlyL  
161 AGCAGGGGGC CCGGATGAGC CTGCTGGGGA AGCCGCTCTC TTACACGAGT AGCCAGAGCT GCCGGCGCAA CGTCAAGTAC  
ysGlnGlyAla ArgMETSer LeuLeuGlyLys ProLeuSer TyrThrSer SerGlnSerCys ArgArgAsn VallysTyr  
241 CGCGGGGTGC AGAACTACCT GTACAACGTG CTGGAGAGAC CCCGCGGCTG GCGGTTCATC TACCACGCTT TCGTTTTTCT  
ArgArgValGln AsnTyrLeu TyrAsnVal LeuGluArgPro ArgGlyTrp AlaPheIle TyrHisAlaPhe ValPheLe  
321 CCTGTCTTTT GGTGCTTGA TTTTGTCACT GTTTCTACC ATCCTGAGC ACACAAATTT GGCCTCAAGT TGCCTCTTGA  
uLeuValPhe GlyCysLeuIle LeuSerVal PheSerThr IleProGluHis ThrLysLeu AlaSerSer CysLeuLeuI  
401 TCCTGGAGTT CGTGATGATT GTCGCTCTTG GTTTGGAGTT CATCATTCGA ATCTGGTCTG CCGGTTGCTG TTGTGCATAT  
leLeuGluPhe ValMETile ValValPheGly LeuGluPhe IleIleArg IleTrpSerAla GlyCysCys CysArgTyr  
481 AGAGGATGGC AAGGAAGACT GAGGTTTGCT CGAAAGCCCT TCTGTGTTAT AGATACCATT GTTCTTATCG CTTCAATAGC  
ArgGlyTrpGln GlyArgLeu ArgPheAla ArgLysProPhe CysValIle AspThrIle ValLeuIleAla SerIleAl  
561 AGTTGTTTCT GCAAAAATC AGGGTAATAT TTTTGCCACG TCTGCACTCA GAAGTCTCCG TTTCTACAG A  
aValValSer AlaLysThrGln GlyAsnIle PheAlaThr SerAlaLeuArg SerLeuArg PheLeuGln

Has open reading frame.

HB\_Q5\_contig\_PCR\_no\_116\_M13rev.  
[Strand]

EcoRV stop

1 GAATTCGGCT TAACTAGAA CTTATTTTCAG TTTGACATGA GGCAAGCTGA GGGCATCTGT ACTTTCTCCT GCCTTACAAA  
81 TGCTCTGAGA TGATCGTGAC CTTCAGTCC TTAGAGAGTC TGATGCAAAG GCAGCTTCCC TGGCAGGCTG CCGTGCGGCA  
25 161 TCAAAAGTGT CTGTCTCTGT CTCTCGGGA CCCACCTCTT CATCAGTTAT AAACAATTTG GATTCCCTCC ATTTGGGGTA  
241 AAAATCTTGG CTGCCTCTGG AGCCACTTGA CTCACTCCCT GAAAGTTGTA TATTCAGTTC CTCGGTTCGAC CTGATCAGGT  
321 TTTGCACAGA CAAAGATTG CCCAAGTCCT TCGGCACCAT GGGACAGACA GACAACAGAG TTTCTCTCTCC CATGTCAAAG  
401 CTTTCTCTCA TAGAACGGTC CTTGGTGAGA TTTGACTGTG CAACCTGACA TTT

#116, in TA - in wrong orientation but can cut out with  
EcoRV at the 5' end of insert, and XbaI at the 3' end of  
the insert. These sites were designed with the primer used  
in PCR:

35

SIGNED

DATE

WITNESSED AND  
UNDERSTOOD BY:

DATE

CROSS REFERENCES:

Cont'd →

DATE 1/25/85 PROJ. NO. \_\_\_\_\_ EXPT. NO. \_\_\_\_\_  
 SUBJECT hKEN25

→ cont'd

INTA

HB\_Q5\_contig\_PCR\_no\_116\_M13rev.  
 [Strand]

rev'd / comp'd


1 AAATGTCAGG TTGCACAGTC AAATCTCACC AAGGACCGTT CTATGAGGAA AAGCTTTGAC ATGGGAGGAG AAACCTCTGTT  
 LysCysGlnVal AlaGlnSer AsnLeuThr LysAspArgSer METArgLys SerPheAsp METGlyGlyGlu ThrLeuLe  
 81 GTCTGTCTGT CCCATGGTGC CGAAGGACTT GGGCAAATCT TTGTCTGTGC AAAACCTGAT CAGGTCGACC GAGGAACCTGA  
 uSerValCys ProMETValPro LysAspLeu GlyLysSer LeuSerValGln AsnLeuIle ArgSerThr GluGluLeuA  
 161 ATATACAACCT TTCAGGGAGT GAGTCAAGTG GCTCCAGAGG CAGCCAAGAT TTTTACCCCA AATGGAGGGA ATCCAAATTG  
 snIleGlnLeu SerGlySer GluSerSerGly SerArgGly SerGlnAsp PheTyrProLys TrpArgGlu SerLysLeu  
 241 TTTATAACTG ATGAAGAGGT GGGTCCCGAA GAGACAGAGA CAGACACTTT TGATGCCGCA CCGCAGCCTG CCAGGGAAGC  
 PheIleThrAsp GluGluVal GlyProGlu GluThrGluThr AspThrPhe AspAlaAla ProGlnProAla ArgGluAl  
 321 TGCCTTTGCA TCAGACTCTC TAAGGACTGG AAGGTCACGA TCATCTCAGA GCATTTGTAA GGCAGGAGAA AGTACAGATG  
 aAlaPheAla SerAspSerLeu ArgThrGly ArgSerArg SerSerGlnSer IleCysLys AlaGlyGlu SerThrAspA  
 401 CCCTCAGCTT GCCTCATGTC AAACCTGAAAT AAGTCTCTAGA TTAAGCCGAA TTC  
 laLeuSerLeu ProHisVal LysLeuLys STP ValLeuAsp STP AlaGlu Phe

Has open reading frame up to the stop site.

Digest ~5 µg of clone #116 (H<sub>B</sub> Brn):

DNA NEB 2 10X BFA H<sub>2</sub>O Enz

20 10X 5 5 26 2 EcoRV, 2 XbaI, 37°, 2 hrs. Run on gel:

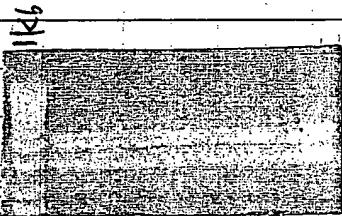

~2.7 kb  
 (~2 µg of insert)

Cut out indicated band, clean up w/  
 Clinch's Nucleotrap kit. Elute w/ H<sub>2</sub>O  
 (50%) - use 3λ to ligate into  
 pCDNA3.1(+).Neo (5' EcoRV / 3' XbaI) -  
 16° o.n.

3/9: Transform 2λ into DH5 & competent cells. Plate out 25λ, 100λ into  
 100 µg/ml Amp plates, 37° o.n. Setup miniprep:

Results: Self-lig (25λ) ~ 100 colonies

Expt'l (25λ) ~ 250 colonies #151-160

Wizard Plus. Digest 2λ w/ EcoRV/XbaI, run on gel:

OD<sub>260</sub> 1.50:

# 260.0

S# For sequencing results see  
 Book 47451, page 24

155 0.1532 380 ng/λ  
 156 0.1658 415 ng/λ  
 157 0.1356 340 ng/λ

Cont'd →

SIGNED

DATE

WITNESSED AND  
 UNDERSTOOD BY:

DATE

DATE: \_\_\_\_\_

PROJ. NO. \_\_\_\_\_

EXPT. NO. \_\_\_\_\_

SUBJECT

hcnos

→ antd

(See Q5 sequence on pp. 186-188 for details)

Design primers for sequencing entire Q5.

→

primers

highlighted  
in sequenceHHMI Biopolymer/Keck Foundation Biotechnology Resource Laboratory  
Yale University School of Medicine

## Oligonucleotide Synthesis Facility

## OLIGO ORDER DETAILS

DO NOT WRITE IN THIS MARGIN

1 Order Date: Q5/396f  
 Sequence Name: 5690  
 Length: 22mer  
 Scale: 40.0  
 Sequence: 5'-GGAGTTCGTGATGATTGTCGTC-3'  
 Notes:  
 Synthesis Date: Chargaff

Order Date: Q5/875f  
 Sequence Name: 5691  
 Length: 22mer  
 Scale: 40.0  
 Sequence: 5'-GGTGGGGCACAATTACATTGAC-3'  
 Notes:  
 Synthesis Date: Chargaff

2 Order Date: Q5/1242f  
 Sequence Name: 5692  
 Length: 19mer  
 Scale: 40.0  
 Sequence: 5'-GACAAGCCTCAGTAGGTGA-3'  
 Notes:  
 Synthesis Date: Chargaff

3 Order Date: Q5/1626f  
 Sequence Name: 5693  
 Length: 21mer  
 Scale: 40.0  
 Sequence: 5'-GCCTTCAAACACGTGTTGATC-3'  
 Notes:  
 Synthesis Date: Chargaff

Order Date: Q5/2078f  
 Sequence Name: 5694  
 Length: 23mer  
 Scale: 40.0  
 Sequence: 5'-CCAATTAGTCAAAGCGATGGCTC-3'  
 Notes:  
 Synthesis Date: Chargaff

Order Date: Q5/2478f  
 Sequence Name: 5695  
 Length: 21mer  
 Scale: 40.0  
 Sequence: 5'-AGTCAAGTGGCTCCAGAGGCA-3'  
 Notes:  
 Synthesis Date: Chargaff

35

Continued in Book 47451, page 21

SIGNED

DATE

WITNESSED AND  
UNDERSTOOD BY:

DATE

CROSS REFERENCES: